

**PRODUCT CHARACTERISTICS**


<b>Concentration:</b>	1-1.5 mg/mL determined by UV absorbance at 280 nm (see vial label)
<b>Specificity:</b>	Reacts to HIV-1 p24 protein
<b>Source:</b>	Murine monoclonal IgG1 from serum-free culture supernatant
<b>Purification:</b>	Protein G chromatography
<b>Purity:</b>	≥ 90% (Coomassie blue SDS-PAGE)
<b>Immunogen:</b>	HIV-1 (Strain IIIB) lysate
<b>Formulation:</b>	1 mg in PBS, pH 7.2 No preservatives

**STORAGE**

Store at -10 °C or below. Repetitive freezing and thawing are not recommended (aliquot as necessary). Thawed material may be stored at 4°C for short-term usage.

**APPLICATIONS**

**Western Blot:** Use at ~1:1000 dilution

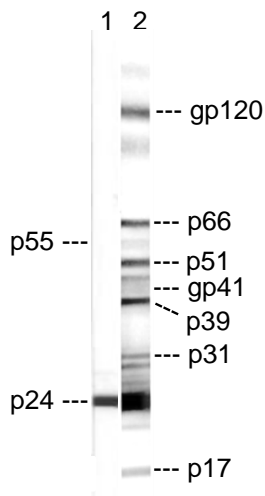

HIV-1 IIIB virus lysate proteins were separated by electrophoresis and transferred onto a nitrocellulose membrane. The following antibodies were used in Western blot analysis.

Strip 1 - Anti-HIV-1 p24 Clone 38/8.7.47

Strip 2 - Human Anti-HIV polyclonal

Additional proteins bands observed on strip 1 are derived from the precursor HIV-1 gag p55 and degradation products.

**ELISA:** Use at ~1:1000 dilution

Recommended dilutions are intended as a guide. Optimal dilutions should be determined by each laboratory based on the application. Conditions for applications such as immunoprecipitation, immunofluorescence and other EIA assays must be determined experimentally by the investigator.

**For Research Use Only NOT for in vitro Diagnostic Use**

PI0801125  
 Revision: 04  
 Effective Date: 02/21/2024

<b>REF</b>	Catalog Number		Temperature Limitation
<b>LOT</b>	Batch Code		Expiration Date
<b>RUO</b>	For Research Use Only		Biological Risk
	Manufacturer		

PCA# 24-059  
Page 1 of 1