

PRODUCT CHARACTERISTICS


Concentration:	1-1.5 mg/mL determined by UV absorbance at 280 nm (see vial label)
Specificity:	Reacts to HIV-1 p24 protein, may react to some strains of SIV and HIV-2
Source:	Murine monoclonal IgG1 from serum-free culture supernatant
Purification:	Protein G chromatography
Purity:	≥ 90% (Coomassie blue SDS-PAGE)
Immunogen:	HIV-1 (Strain IIIB) lysate
Formulation:	1 mg in PBS, pH 7.2 No preservatives

STORAGE

Store at -10 °C or below. Repetitive freezing and thawing are not recommended (aliquot as necessary). Thawed material may be stored at 4°C for short-term usage.

APPLICATIONS

Western Blot: Use at ~1:1000 dilution

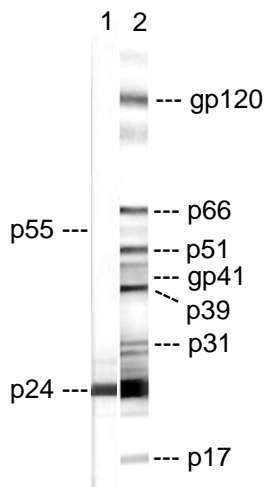

HIV-1 IIIB virus lysate proteins were separated by electrophoresis and transferred onto a nitrocellulose membrane. The following antibodies were used in Western blot analysis.

Strip 1 - Anti-HIV-1 p24 Clone 39/6.14

Strip 2 - Human Anti-HIV polyclonal

Additional proteins bands observed on strip 1 are derived from the precursor HIV-1 gag p55 and degradation products.

ELISA: Use at ~1:1000 dilution

Recommended dilutions are intended as a guide. Optimal dilutions should be determined by each laboratory based on the application. Conditions for applications such as immunoprecipitation, immunofluorescence and other EIA assays must be determined experimentally by the investigator.

For Research Use Only NOT for in vitro Diagnostic Use

PI0802014
Revision: 05
Effective Date: 02/21/2024

REF	Catalog Number		Temperature Limitation
LOT	Batch Code		Expiration Date
RUO	For Research Use Only		Biological Risk
	Manufacturer		

PCA# 24-059
Page 1 of 1